User's Manual

Protector[®] Stainless Steel Radioisotope Laboratory Fume Hoods

Models

1204 Series

1205 Series

1206 Series

1208 Series

Protector

To receive important product updates,
complete your product registration card
online at register.labconco.com

*Protecting your
laboratory environment*

LABCONCO[®]

Labconco Corporation
8811 Prospect Avenue
Kansas City, MO 64132-2696
800-821-5525, 816-333-8811
FAX 816-363-0130
E-MAIL labconco@labconco.com
HOME PAGE www.labconco.com

Please read the User's Manual before operating the equipment.

The information contained in this manual and the accompanying products are copyrighted and all rights reserved by Labconco Corporation. Labconco Corporation reserves the right to make periodic design changes without obligation to notify any person or entity of such change.

Warranty

Labconco provides a warranty on all parts and factory workmanship. The warranty includes areas of defective material and workmanship, provided such defect results from normal and proper use of the equipment.

The warranty for all Labconco products will expire one year from date of installation or two years from date of shipment from Labconco, whichever is sooner, except the following;

- Purifier® Logic®+ Biological Safety Cabinets and PuriCare® Procedure Stations carry a five-year warranty from date of installation or six years from date of shipment from Labconco, whichever is sooner.
- SteamScrubber® & FlaskScrubber® Glassware Washers carry a two-year warranty from date of installation or three years from date of shipment from Labconco, whichever is sooner.
- Blood Drawing Chairs carry a ten year warranty.
- Carts carry a lifetime warranty.
- Glassware is not warranted from breakage when dropped or mishandled.

This limited warranty covers parts and labor, but not transportation and insurance charges. In the event of a warranty claim, contact Labconco Corporation or the dealer who sold you the product. If the cause is determined to be a manufacturing fault, the dealer or Labconco Corporation will repair or replace all defective parts to restore the unit to operation. Under no circumstances shall Labconco Corporation be liable for indirect, consequential, or special damages of any kind. This statement may be altered by a specific published amendment. No individual has authorization to alter the provisions of this warranty policy or its amendments. Lamps and filters are not covered by this warranty. Damage due to corrosion or accidental breakage is not covered.

Returned or Damaged Goods

Do not return goods without the prior authorization from Labconco. Unauthorized returns will not be accepted. If your shipment was damaged in transit, you must file a claim directly with the freight carrier. Labconco Corporation and its dealers are not responsible for shipping damages.

The United States Interstate Commerce Commission rules require that claims be filed with the delivery carrier within fifteen (15) days of delivery.

Limitation of Liability

The disposal and/or emission of substances used in connection with this equipment may be governed by various federal, state, or local regulations. All users of this equipment are required to become familiar with any regulations that apply in the user's area concerning the dumping of waste materials in or upon water, land, or air and to comply with such regulations. Labconco Corporation is held harmless with respect to user's compliance with such regulations.

Contacting Labconco Corporation

If you have questions that are not addressed in this manual, or if you need technical assistance, contact Labconco's Customer Service Department or Labconco's Product Service Department at 1-800-821-5525 or 1-816-333-8811, between the hours of 7:00 a.m. and 6:00 p.m., Central Standard Time.

TABLE OF CONTENTS

CHAPTER 1: INTRODUCTION	1
About This Manual	2
Typographical Conventions (Conventions Typographiques)	3
Your Next Step	5
CHAPTER 2: PREREQUISITES	6
Location Requirements	7
Support Requirements	7
Exhaust Requirements	7
Filter Pack Requirements	8
Electrical Requirements (Exigences électriques)	8
Service Line Requirements	9
Space Requirements	9
CHAPTER 3: GETTING STARTED	10
Unpacking Your Fume Hood	11
Removing the Shipping Skid	11
Sash Weight Release	12
Installing the Hood on a Supporting Structure & Hardboard Work Top	12
Connecting to the Hood Exhaust System	15
Connection to the Optional Filter Packs	16
Connecting the Electrical Supply Source	16
Connecting the Service Lines to the Protector Hood	19
Sealing the Protector Hood to the Work Surface	20
Certifying the Protector Radioisotope Fume Hood	20
CHAPTER 4: PERFORMANCE FEATURES AND SAFETY PRECAUTIONS	21
Performance Features	21
Safety Precautions (Consignes de sécurité)	25

CHAPTER 5: USING THE PROTECTOR RADIOISOTOPE	28
Operating the Vertical-Rising Glass Sash	28
Operating the Blower	28
Operating the Lights	28
Operating the Washdown Control Valve	28
Working in Your Protector Radioisotope Fume Hood	29
 CHAPTER 6: MAINTAINING THE PROTECTOR RADIOISOTOPE	 31
Service Safety Precautions (Mesures de sécurité des services)	31
Routine Maintenance Schedule	32
Routine Service Operations	33
 CHAPTER 7: MODIFYING THE PROTECTOR RADIOISOTOPE	 34
Installing Hardboard Support Surface	34
Installing Ceiling Enclosures Above the Fume Hood	35
Installing Rear Panels Behind the Fume Hood	35
Installing Guardian Airflow Monitors	35
Sash Stop Kit – Field Installation	35
Installing an Electrical Duplex Outlet (Installation d'une sortie Duplex électrique)	36
 CHAPTER 8: TROUBLESHOOTING	 37
 APPENDIX A: PROTECTOR RADIOISOTOPE COMPONENTS	 40
 APPENDIX B: PROTECTOR RADIOISOTOPE DIMENSIONS	 43
 APPENDIX C: PROTECTOR RADIOISOTOPE SPECIFICATIONS (ANNEXE C: PROTECTEUR DE RADIO-ISOTOPES SPECIFICATIONS)	 47
 APPENDIX D: SERIAL NUMBER TAG DESCRIPTION (ANNEXE D: DESCRIPTION DE SÉRIE TAG NUMBER)	 50
 APPENDIX E: REFERENCES	 56
 DECLARATION OF CONFORMITY	 59

CHAPTER 1

INTRODUCTION

Congratulations on your purchase of a Labconco Protector® Stainless Steel Radioisotope Laboratory Fume Hood. The Protector Stainless Steel Radioisotope Laboratory Hood is designed to efficiently ventilate procedures using radiochemicals. It is the result of Labconco's more than 50 years experience in manufacturing fume hoods, and users like you suggested many of its features to us.

The Labconco Protector Stainless Steel Radioisotope Fume Hood has been engineered to provide maximum safety in a laboratory while using low level Alpha and lower level Beta particles, and effectively contain toxic, noxious, or other harmful materials when properly installed. The Protector Radioisotope Hood offers many unique features to enhance safety, performance, and visibility. To take full advantage of them, please acquaint yourself with this manual and keep it handy for future reference. If you are unfamiliar with how fume hoods operate, please review *Chapter 4: Performance Features and Safety Precautions* before you begin working in the fume hood. Even if you are an experienced fume hood user, please review *Chapter 5: Using the Protector Radioisotope*, which describes your Protector Radioisotope Hood's features so that you can use the hood efficiently.

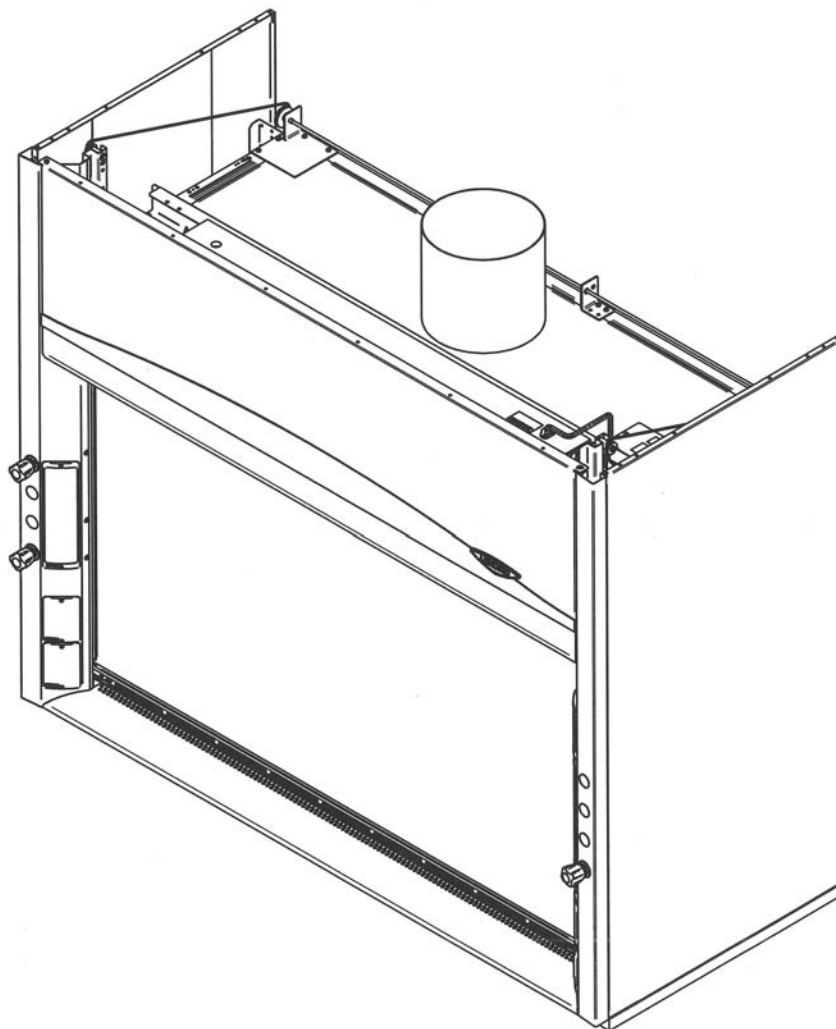
Figure 1-1

About This Manual

This manual is designed to help you learn how to install, use, and maintain your laboratory fume hood. Instructions for installing optional equipment on your hood are also included.

Chapter 1: Introduction provides a brief overview of the laboratory fume hood, explains the organization of the manual, and defines the typographical conventions used in the manual.

Chapter 2: Prerequisites explains what you need to do to prepare your site before you install your laboratory fume hood. Electrical and service requirements are discussed.

Chapter 3: Getting Started contains the information you need to properly unpack, inspect, install, and certify your laboratory fume hood.

Chapter 4: Performance Features and Safety Precautions explains how the Protector operates and the appropriate precautions you should take when using the fume hood.

Chapter 5: Using The Protector Radioisotope discusses the basic operation of your fume hood. Information on how to prepare, use and shut down your Protector Hood are included.

Chapter 6: Maintaining The Protector Radioisotope explains how to perform routine maintenance on your fume hood.

Chapter 7: Modifying The Protector Radioisotope explains how to modify the fume hood or add accessories.

Chapter 8: Troubleshooting contains a table of problems you may encounter while using your laboratory fume hood including the probable causes of the problems and suggested corrective actions.

Appendix A: Protector Radioisotope Components contains labeled diagrams of all of the components of the fume hoods.

Appendix B: Protector Radioisotope Dimensions contains comprehensive diagrams showing all of the dimensions for the laboratory fume hoods.

Appendix C: Protector Radioisotope Specifications contains the electrical requirements for laboratory fume hood. Wiring diagrams are also included.

Appendix D: Serial Number Tag Description provides current rating code used on serial number tag.

Appendix E: References lists the various resources available that deal with laboratory fume hoods.

Typographical Conventions

Recognizing the following typographical conventions will help you understand and use this manual:

- Book, chapter, and section titles are shown in italic type (e.g., *Chapter 3: Getting Started*).
- Steps required to perform a task are presented in a numbered format.
- Comments located in the margins provide suggestions, reminders, and references.

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
- Critical information is presented in boldface type in paragraphs that are preceded by the wrench icon. These operations should only be performed by a trained certifier or contractor. Failure to comply with the information following a wrench icon may result in injury to the user or permanent damage to your hood.

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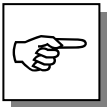- The S icon indicates the text is specific to the standard model.

- The A icon indicates the text is specific to the A-Style Combination Sash Model.

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- CAUTION – Hot Surface.
- AVERTIR – Surface Chaude

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Your Next Step

If your Fume Hood needs to be installed, proceed to *Chapter 2: Prerequisites* to ensure your installation site meets all of the requirements. Then, go to *Chapter 3: Getting Started* for instructions on how to install your laboratory fume hood and make all of the necessary connections.

If you would like to review how laboratory fume hoods operate, go to *Chapter 4: Performance Features and Safety Precautions*.

For information on the operational characteristics of your laboratory fume hood, go to *Chapter 5: Using The Protector Radioisotope*.

If your laboratory fume hood is installed and you need to perform routine maintenance on the cabinet, proceed to *Chapter 6: Maintaining Your Protector Radioisotope*.

For information on making modifications to the configuration of your fume hood, go to *Chapter 7: Modifying The Protector Radioisotope*.

Refer to *Chapter 8: Troubleshooting* if you are experiencing problems with your fume hood.

CHAPTER 2

PREREQUISITES

Before you install your laboratory fume hood, you need to prepare your site for installation. Carefully examine the location where you intend to install your hood. You must be certain that the area is level and of solid construction. In addition, a dedicated source of electrical power must be located near the installation site.

Carefully read this chapter to learn the requirements for your installation site:

- The location requirements.
- The support requirements.
- The exhaust requirements.
- The optional filter pack requirements.
- The electrical power requirements.
- The service line requirements.
- The space requirements.

Refer to *Appendix B: Protector Radioisotope Dimensions* for complete fume hood dimensions.

Refer to *Appendix C: Protector Radioisotope Specifications* for complete laboratory fume hood electrical and environmental conditions, specifications and requirements.

Location Requirements

The fume hood should be located away from traffic patterns, doors, windows, fans, ventilation registers, and any other air-handling device that could disrupt its airflow patterns. All windows in the room should be closed.

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Support Requirements

DO NOT install the fume hood on a cart, dolly, or mobile bench. ALL Protector Hood installations must be permanent and stationary. The supporting structure usually consists of a base cabinet and chemically resistant work surface.

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Exhaust Requirements

The exhaust duct connection has been designed for 11" nominal duct (11.00" OD) to allow for minimum static pressure loss while operating at 100 fpm face velocities. The 11" diameter exhaust duct also allows for proper transport velocities away from the hood in the 1000 fpm to 2500 fpm range. The exhaust duct used in conjunction with the Radioisotope Fume Hood should be fabricated out of Type I, unplasticized PVC or Type 304 stainless steel. The proper exhaust volume and static pressure loss are listed next for each hood model:

Chapter 2: Prerequisites

Airflows and Static Pressure

Hood Type	Eco Energy Saving Air Foil					
	Face	Velocity (fpm)	Airflow Volumetric Rate (CFM) @ Static Pressure (inches of water)			
	Sash at Full open at 28"	Sash at 62.5% open at 18"	4' Hood	5' Hood	6' Hood	8' Hood
Stainless Hoods						
Perchloric and Radioisotope	125	200	905, 0.39"	1195, 0.52"	1475, 0.78"	2050, 0.44"
	100	160	725, 0.25"	955, 0.33"	1180, 0.50"	1640, 0.28"
	N/A	125	565, 0.15"	745, 0.20"	920, 0.31"	1280, 0.17"
11" OD Exhaust	N/A	100	455, 0.10"	595, 0.13"	740, 0.20"	1025, 0.11"

Recommendations highlighted

Proper Radioisotope blower selection can be determined from these exhaust requirements and the total system static pressure loss. Contact Labconco Customer Service for assistance in sizing a blower system.

Filter Pack Requirements

Contact your local EPA or NRC office to determine if exhaust filtration is required. If required, these filter packs should be located as close as possible to the hood to expose as little of the ductwork system as possible to the radiochemicals being filtered. Consult your Safety Officer directly for specific filter ratings on radiochemicals being used.

Electrical Requirements

The Protector Hood models feature internal wiring for the fluorescent light assembly and light switch. All internal wiring is terminated at the single point wiring junction box for hook-up by a qualified electrician. The blower switch, and light switch wires are also terminated at the single point wiring junction box for hook-up by a qualified electrician. Refer to *Chapter 3: Getting Started* and *Appendix C: Protector Radioisotope Specifications* for the wiring diagram for proper electrical installation.

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Service Line Requirements

All service lines to the laboratory fume hood should be ¼ inch outside diameter, copper (brass for natural gas), and equipped with an easily accessible shut-off valve, should disconnection be required. If the service line pressure exceeds 40 PSI, it must be equipped with a pressure regulator to reduce the line pressure. Please check with local codes for other requirements.

Space Requirements

The dimensions for the different models are shown in *Appendix B: Protector Radioisotope Dimensions*.

CHAPTER 3

GETTING STARTED

Now that the site for your laboratory fume hood is properly prepared, you are ready to unpack, inspect, install, and certify your unit. Read this chapter to learn how to:

- Unpack and move your Protector Radioisotope Hood.
- Set up the fume hood with the supporting structure and hardboard work top.
- Connect to an exhaust system suitable for radiochemicals.
- Connect the optional filter packs to the exhaust system (not sold through Labconco).
- Connect the electrical supply source.
- Connect the service lines.
- Arrange certification of your Protector Radioisotope Hood.

Depending upon which model you are installing, you may need common plumbing and electrical installation tools in addition to 5/16", 3/8", 7/16", and 1/2" wrenches, ratchets, sockets, a nut driver set, a flat-blade screwdriver, a Phillips screwdriver, and a carpenter level to complete the instructions in the chapter.


The Protector Radioisotope Hood models weigh between 400 to 800 lbs. (182-363 kg). The shipping skid allows for lifting with a mechanical lift truck or floor jack. If you must lift the fume hood manually, follow safe-lifting guidelines. Normally, the fume hood can be slid off a hydraulic lift table and be placed into position on top of the work surface. Do not lift by the front air foil.

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
Unpacking Your Fume Hood

Carefully remove the shrink-wrap or carton on your fume hood and inspect it for damage that may have occurred in transit. If your unit is damaged, notify the delivery carrier immediately and retain the entire shipment intact for inspection by the carrier.

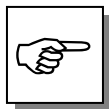DO NOT RETURN GOODS WITHOUT THE PRIOR AUTHORIZATION OF LABCONCO. UNAUTHORIZED RETURNS WILL NOT BE ACCEPTED.

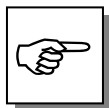IF YOUR HOOD WAS DAMAGED IN TRANSIT, YOU MUST FILE A CLAIM DIRECTLY WITH THE FREIGHT CARRIER. LABCONCO CORPORATION AND ITS DEALERS ARE NOT RESPONSIBLE FOR SHIPPING DAMAGES.

The United States Interstate Commerce Commission rules require that claims be filed with the delivery carrier within fifteen (15) days of delivery.


Do not discard the shipping skid or packing material for your fume hood until you have checked all of the components and installed and tested the unit. **The Perchloric Acid fume hood baffle is installed as shown in Figure 1-1.** Do not remove the fume hood from its shipping skid until it is ready to be placed into its final location. Move the unit by placing a flat, low dolly under the shipping skid, or by using a floor jack.

Do not move the hood by tilting it onto a hand truck.

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Removing the Shipping Skid

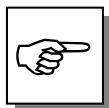LEAVE THE FUME HOOD ATTACHED TO ITS SHIPPING SKID UNTIL IT IS AS CLOSE TO ITS FINAL LOCATION AS POSSIBLE. MOVE THE HOOD BY USING A SUITABLE FLOOR JACK, OR BY PLACING A FURNITURE DOLLY UNDERNEATH THE SKID. DO NOT MOVE THE HOOD BY TILTING IT ONTO A HAND TRUCK.

After you verify the fume hood components, move your hood to the location where you want to install it. Follow the steps listed next to remove the shipping skid from your unit.


1. Remove the side panels by unscrewing the Phillips screws.

Chapter 3: Getting Started

2. Find the hardware (bolts, washers, nuts) that attach the fume hood to the skid and remove the hardware. Some hardware is on the sides and some is on the back.

Sash Weight Release

To protect the fume hood from damage in shipment, the sash weight has been secured to the back of the fume hood with screws. Simply remove the screws and make sure the sash cables or chains are on the pulleys or sprockets before operation of the sash.

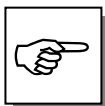NOTE: THE SASH WEIGHT ITSELF WAS INDIVIDUALLY MATCHED FOR THIS SPECIFIC HOOD AND SHOULD NOT BE EXCHANGED ON ANY OTHER UNIT.

Installing the Hood on a Supporting Structure and Hardboard Work Top

The Protector Hood is heavy! Use caution when lifting or moving the unit.

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When installing the Protector Radioisotope Hood onto the lower hardboard work top, ensure that the structure can safely support the combined weight of the fume hood and any related equipment. **The hardboard work top is aligned flush with the back of the fume hood for good airflow: this will provide the correct spacing under the air foil for proper bypass airflow. The lower base cabinets are placed flush with the front of the hardboard work top as shown in Figure 3-2. The hardboard work top supports the Radioisotope Hood.**

WARNING: It is important to support the rear of the work surface and fume hood. The cross support provides support for the bottom of the work surface. Install the cross support after the base cabinets and work surface are leveled and before installing the hood.

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The following are instructions for mounting a cross support:

1. Level the base cabinets and the hardboard work top. Work top should be placed flush with the back of the fume hood as shown in Figure 3-1.
2. Scribe a line on the wall or back of the base cabinet to locate the support under the work surface.
3. Mount the support by attaching it to the wall or base cabinet.
4. Anchor cabinet(s) with seismic restraints if required by the local code in your area.

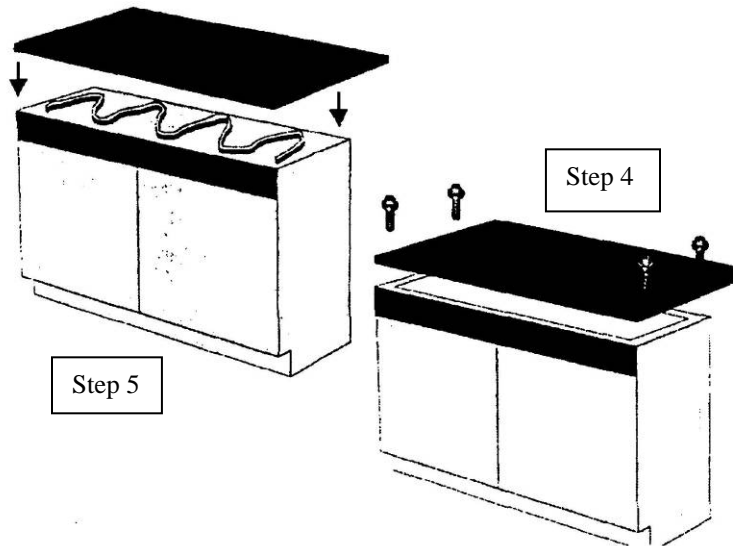
Figure 3-1

5. Secure the hardboard surface to the top of the base cabinet(s) using adhesives with the Storage Cabinet(s). Be sure the Radioisotope Hood drain trough is accounted for in the relief area of the hardboard top. See figure 3-2.
6. After completing step 5, apply a sufficient quantity of adhesive to the top surface of the hardboard.
7. Place the hood on top of the work surface and cross support.
8. Place heavy weights on the inside of the fume hood for approximately 24 hours. This provides a proper bond of the adhesive. Protect the hood's integral work surface from possible damage by wrapping any weighted material first.
9. Place the hood on top of the work top and cross support.

Connecting to the Hood Exhaust System

WARNING: The weight of the exhaust ductwork system must be supported independently of the hood superstructure. Do not allow this weight to be supported by the hood structure as damage to the hood may occur.

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The exhaust connection should be installed by a qualified HVAC contractor. The exhaust connection on your hood is manufactured with 11" Stainless Steel Type 304 nominal pipe (11.00" OD) to allow for minimum static pressure loss with proper transport velocities away from the hood. Consult Labconco Customer Service should you require help sizing your blower for the exhaust volume and total system static pressure loss.

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The selected exhaust duct material must be fabricated out of type 1 unplasticized polyvinyl chloride or type 304 stainless steel material. The duct should be fully welded in the case of the stainless steel material or solvent bonded in the case of polyvinyl chloride so that it does not have any cracks or holes to allow for the collection of radiochemicals.

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Connecting to the Optional Filter Packs

Both HEPA (High Efficiency Particulate Air) and activated carbon filter pack assemblies are available for use with the Protector Radioisotope Hood. Please contact your local EPA or NRC office for manufacturers of exhaust filters appropriate for use with radio-chemicals. The HEPA filter assembly can be used to remove most particulate forms of radio-chemicals being used today. The activated carbon filter is used in the removal of gaseous form radiochemicals and can be used in conjunction with the HEPA filter, depending on the materials being used inside the fume hood. The filters themselves should be located as close as possible to the hood to expose as little of the ductwork system as possible to the radiochemicals being used. Consult your company Safety Officer directly for specific filter ratings on radiochemicals being used.

IMPORTANT WARNING: Used filters contain hazardous material. Observe all federal, state and local EPA regulations for their handling and disposal.

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Connecting the Electrical Supply Source to the Protector Fume Hood

Prior to connecting any electrical wiring to the fume hood structure, refer to the hood identification plate for the proper electrical requirements of your specific model.

WARNING: The building electrical supply system for Protector Hoods is required to include overload protection. A switch or circuit breaker shall be in close proximity to the equipment and within easy reach of the operator. The switch or circuit breaker is to be marked as the disconnecting device for the equipment.

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The identification plate, model number, serial number, and electrical connection boxes are accessible from the front of the fume hood by removing the front panel.

The Protector Radioisotope Hood is normally wired for 115 Volt, 50/60 Hz, 20 Amp or 230Volt, 50/60 Hz, 10 Amp electrical service. Check the I.D. plate behind the front panel for voltage verification. The number of circuits varies depending on the model. All of the electrical connections are terminated at the field wiring terminal box for hook-up by a qualified electrician. We recommend each circuit be a dedicated branch circuit. However, if wired together the maximum load allowed is the sum of individual outlets plus the rating of the unit (i.e. 2 Amps remote, 8 Amps integral blower) not to exceed the 10A or 20A electrical service. The single point internal junction box is used for the connection of the lights, blower, and duplex outlets. Refer to the wiring diagram for your Protector Radioisotope in *Appendix C: Protector Radioisotope Specifications*.

The fume hood is required to be grounded to the MAINS protective earthing ground for safe operation. Using a ring terminal sized for a 10-24 machine screw, connect the MAINS ground conductor to the grounding lug marked with the protective earthing symbol, . Only MAINS ground conductors should be connected to the protective earthing ground lug, no other conductors should be connected to this grounding lug. Using wire nuts connect the MAINS supply conductors to the fume hood supply wires. Insure that the wires are connected as per the appropriate wire color codes for the input voltage. For 115V Phase (Hot) is black and Neutral is white, for 230V Phase1 is brown and Phase2 is blue. Refer to the wiring diagram for your Protector Radioisotope in *Appendix C: Protector Radioisotope Specifications*.

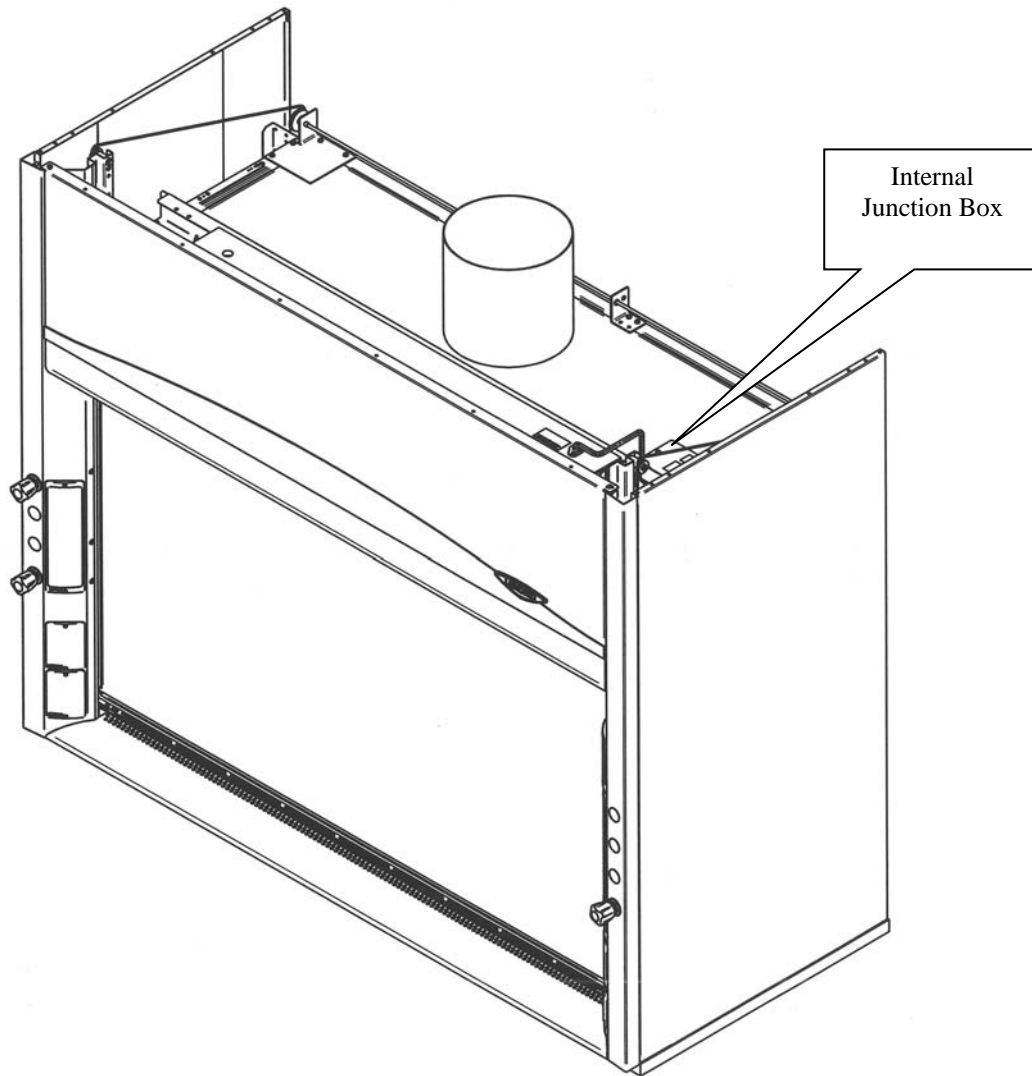
Figure 3-3

All wiring for the fume hood **SHOULD** be performed by a licensed electrician and conform to all local codes. In most cases, the hood will require the use of shielded conduit to protect the wiring into the hood. The grounding connection shall not be made to the terminal box cover.

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The fluorescent light has been mounted outside the top liner panel and is sealed from vapors inside the hood structure. To change the fluorescent light bulbs in your hood, you must first remove the front panel from the hood. Next remove the knock out plugs holding the light fixture in place. Lift fixture up and replace any defective bulbs. Reverse order to reassemble.

Connecting the Service Lines to the Protector Fume Hood

The hoods with service fixtures have been plumbed from the valve to the hose connector or gooseneck for your installation convenience. Supply tubing shall be provided by the qualified installer. Tubing can enter the hood from above, through the back, or through the work surface to make these connections to the service fixtures.

NOTE: Inspect all fittings for leakage. Tighten the fittings slightly if needed.

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CAUTION: Do not use oxygen with any standard service fixture. Contact Labconco Customer Service for oxygen fixture information.

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Should access to the hood plumbing fixture bodies be required, remove the service access plate on the hood front corner posts by loosening their individual screws (see item 11, Figure A-1 in *Appendix A*). The valve body will now be fully exposed for any service work that may be necessary. The service fixtures supplied on your laboratory hood are designed for use with the following services:

- Air
- Cold Water
- Hot Water
- Natural Gas – See caution below
- Vacuum

WARNING: Contact Labconco Customer Service directly before using any service other than those listed above in these valves to assure full compatibility.

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CAUTION: Natural gas should be used only in the service fixture that has been pre-plumbed with brass tubing. Sulfur content of the gas could cause deterioration of standard copper supply lines.

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Sealing the Protector Hood to the Work Surface

When the hood has been set in place, ducted, wired, and plumbed, it should be sealed to the hardboard work top. Materials such as silicone sealants are recommended to seal the hood structure.

Certifying the Protector Radioisotope Fume Hood

The combination of your laboratory hood, exhaust ductwork, and exhaust blower gives you the flexibility to change the airflow at the sash opening of your hood. To determine the actual face velocity at the sash opening, airflow velocity readings will need to be taken. This should be done across the sash opening of the hood in accordance with the *Industrial Ventilation Manual* section on laboratory hoods (see *Appendix E: References*). Labconco recommends an average face velocity at the sash opening of 100 feet per minute. Consult Chapter 2 for proper airflow volumes for your particular model.

Your Protector Fume Hood has been tested at the factory per ASHRAE 110-1995. All hoods achieve an “as manufactured rating” of less than 0.05 part per million (ppm) at 4 liters per minute (lpm); AM<0.05 (consult Labconco for individual fume hood ratings). For “field use” ASHRAE testing contact Labconco Sales Engineering Team or Customer Service for a certified on-site contractor.

NOTE: Face velocity profiles and smoke testing should be done periodically to ensure safe performance.

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CHAPTER 4

PERFORMANCE FEATURES AND SAFETY PRECAUTIONS

Performance Features

The Protector Radioisotope Laboratory Hood has been engineered to provide you the maximum in work area and personal safety. The by-pass air configuration on the Radioisotope Fume Hood is totally dependent on air being supplied to the hood from its surrounding environment for proper operation.

The hood itself is designed to handle radiochemical operations of all types. The hood liner should be decontaminated manually by the operator. Water washdown or flooding the surface with water is standard practice and results in about 50% reduction of contamination on the 304 stainless steel liner material. A secondary means of surface decontamination involves the use of organic solvents that can result in reductions of contamination of up to 75%.

Even though the stainless steel seams are smooth to the touch, there are microscopic voids within the surface in which radioactive particles can lodge and are not easily removed by the above mentioned surface decontamination procedures. Abrasion decontamination, which consists of sanding and/or scouring may be required to dislodge these particles from your hood liner.

Radioisotope materials are, in essence, high-energy particles with different levels of intensity. These levels are, alpha, beta, and gamma. Radioactive particles do not just simply strike a barrier and bounce off, but rather are absorbed by the barrier and physical decomposition on the molecular level can occur. In the case of the alpha and some lower level beta radiation particles, human skin can act as an effective barrier to this radiation without damage to the underlying tissues.


Chapter 4: Performance Features & Safety Precautions

However, in the case of inhalation, low-level emitters can attack inner organ tissues posing a health hazard. Therefore, it is very important in cases even where these low level type emitters are being used that effective sources of ventilation be provided.

The radioisotope hood's stainless steel liner material is intended for use with alpha particles, which are least hazardous, and for most beta type particles, which have a slightly higher mass and consequently feature a higher velocity requiring more shielding than the alpha particles. The gamma emitter particles have no mass but an extremely high velocity. Therefore, they are not easily absorbed and are not normally used within this hood enclosure unless it has been properly lined with lead blocks and a lead glass sash assembly. Consult with your company safety officer for additional safety procedures dealing with your specific materials, before operating this hood.

1. **Unique sash provides maximum visibility of 37.5" high while conserving energy by limiting sash travel to 28"**. Vertical-rising sash may be raised from a closed to 28" operating height. Exhaust volume and blower sizing are based on the 28" height. Optional sash stops are available to limit sash height and reduce energy usage. The Perchloric Acid Hood has a glass sash specifically used for Perchloric Acid.
2. **By-pass airflow design** ensures relatively stable face velocities.
3. **Large usable interior work depth and interior height of 48"** provides ample working space.
4. **Baffle** directs airflow to the rear of the interior to provide efficient airflow. The baffle may be removed for cleaning purposes only.
5. **Exterior access cover plates** are removable for easy access to plumbing valves when access through the sides is not available.
6. **Lift-Away™ front panel** provides easy access to electrical wiring, sash weights, and lighting fixtures.
7. **Energy efficient fluorescent lighting** is located behind a laminated safety glass shield mounted to the top of the hood. The factory-wired lighting is serviceable from outside the hood cavity.
8. **Low-mounted, factory-wired light and blower switches** are ADA compliant.
9. **Clean-Sweep™** air foil allows air to sweep the work surface for maximum containment.
10. **Streamlined corner posts** provide maximum visibility and the flexibility to add services; services are most easily added at the factory, because of the PVC liner.
11. **All hoods are factory-prepared for up to 7 service fixtures**. Additional fixtures are available only as a factory special. Liner hole drilling is required to add another fixture in the field.

12. **One duplex electrical receptacle** is mounted on the right corner post. Receptacles are factory-wired to hood single point junction box. Additional electrical receptacles are available as factory or field installed.
13. **Type 304 Stainless Steel interior liner** is seamless with coved interior corners and integral work surface.
14. **Accessory Guardian™ Digital Airflow Monitor or Guardian Jr. Monitor** continuously monitors face velocity. An audio/visual alarm alerts the user to low airflow conditions. The right corner post is factory prepared to accommodate the Guardian Monitor (sold separately).
15. **Optional Energy-Reducing A-Style Combination Sash Models.** These combination sashes allow the operator to use the hood with sashes either half open vertically or horizontally to conserve energy. Optional sash stops prevent raising the vertical sash above the half-open and fully closed positions unless manually defeated by the operator.
16. **Outside frame of epoxy-coated steel and galvanized steel** is durable and corrosion resistant.
17. **Exhaust connection.** The hood features 11" (11.00" OD pipe) Type 304 Stainless Steel exhaust connections sized to allow for a minimum static pressure loss through the hood structure while providing a good transport velocity through the exhaust system.
18. **Solid hardboard support** provided to support the integral Type 316 Stainless Steel bottom.
19. **Optional Ceiling Enclosure Kits** are available for a decorative facade between the hood and the ceiling.
20. **Optional Sash Stops** provide a means of controlling the operating height of the sash.
21. **Built-in Upper Bypass.** The Radioisotope Hood includes an upper bypass block which proves to be useful when working with heat in digestions.
22.  **CAUTION** – Current rating of receptacle is specified in *Appendix D*.
AVERTIR – Classification des prises de courant est spécifié dans l'Annexe D.
23.  **CAUTION** – See Appendix C and D for complete current rating.
AVERTIR – Voir Annexe C et D pour la classification complète du courant.
24.  **CAUTION** – Flammable Gas.
AVERTIR – Gaz Inflammable.

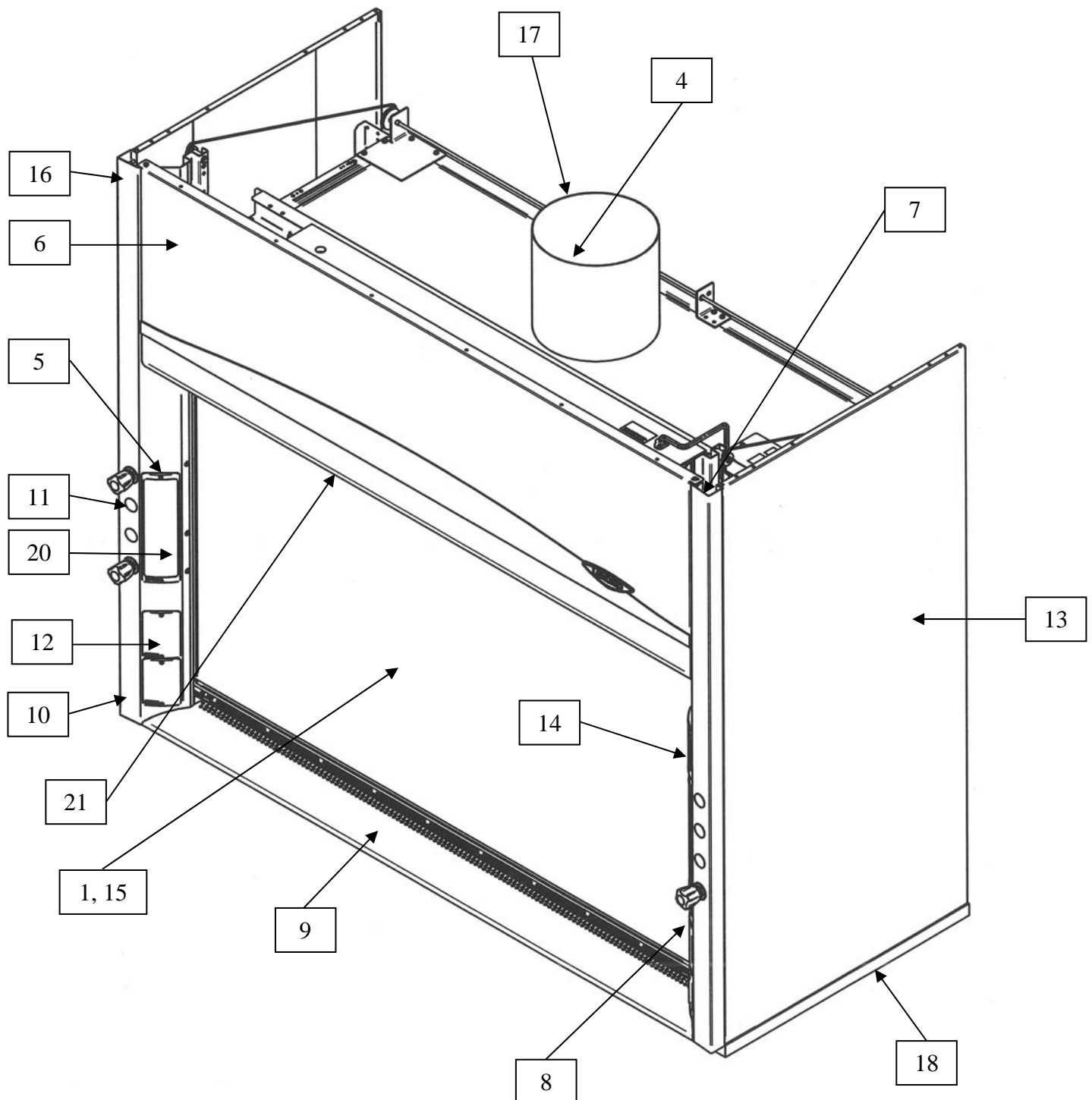
Chapter 4: Performance Features & Safety Precautions**Figure 4-1**

Safety Precautions

Although the laboratory hood has been engineered to maintain optimum operator safety, caution should always be used while working in the hood. Prior to using the hood, check to make sure that the exhaust blower is operating and that air is entering the hood at its specified face velocity.

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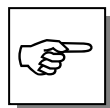USE GOOD HOUSEKEEPING IN THE HOOD AT ALL TIMES. CLEAN UP SPILLS IMMEDIATELY WITH A MILD DETERGENT. PERIODICALLY CLEAN HOOD INTERIOR, INCLUDING FLUORESCENT LIGHT GLASS PANEL. REPLACE BURNED OUT LIGHT BULBS TO MAINTAIN MAXIMUM ILLUMINATION.

DO NOT OVERLOAD THE WORK SURFACE WITH APPARATUS OR WORK MATERIAL. THE SAFE OPERATION OF THE LABORATORY HOOD IS BASED UPON HAVING PROPER AIRFLOW THROUGH THE STRUCTURE. DO NOT PLACE LARGE, BULKY OBJECTS SUCH AS BLOCK HEATERS, DIRECTLY ON THE HOOD WORK SURFACE. INSTEAD, ELEVATE THE OBJECT 2" TO 3" ON BLOCKS TO ALLOW A FLOW OF AIR UNDER THE OBJECT AND INTO THE LOWER REAR BAFFLE EXHAUST SLOT. ENSURE BLOCKS ARE LEVEL AND SECURED IN PLACE.

Blocking the bottom of the baffle at the rear of hood will change the airflow pattern in the hood causing turbulence and possible leakage at the face of the hood. (Don't store containers or supplies against baffles, as this will affect airflow through the hood).

Avoid placing your head inside hood. Keep hands out of hood as much as possible.

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Always work as far back in hood as possible. It is best to keep all chemicals and apparatus 6" inside the front of the hood.

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
Perchloric acid use in Stainless Steel Radioisotope Hoods is not recommended.

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High-level radioisotope materials are not recommended for use in this Radioisotope hood. Only use low-level alpha and lower level beta radiochemical particulates.

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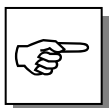AVOID CROSS DRAFTS AND LIMIT TRAFFIC IN FRONT OF THE HOOD. AIR DISTURBANCES CREATED MAY DRAW FUMES OUT OF THE HOOD.

The use of heat-generating equipment in this hood without the exhaust system operating properly can cause damage to the hood.

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The Protector Radioisotope Laboratory Hoods should be certified by a qualified certification technician before initially used. The hood should be re-certified whenever it is relocated, serviced or at least annually thereafter.

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Ensure that the hood is connected to electrical service in accordance with local and national electrical codes. Failure to do so may create a fire or electrical hazard. Do not remove or service any electrical components without first disconnecting the hood from electrical service. Proper operation of the fume hood depends largely upon the hood's location and the operator's work habits. Consult the *Reference Manual in Appendix D*.

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If the unit is not operated as specified in this manual it may impair the protection provided by the unit.

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Do not position the fume hood so that it is difficult to operate the main disconnect device.

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RADIOISOTOPE

S

Operating the Vertical-Rising Glass Sash

Because of the Protector Radioisotope Hood's counterbalanced sash mechanism, it will take only a few pounds of force to move the sash up or down, and you can operate the sash smoothly with one or two hands positioned any where along the handle. The glass sash is suitable for use with Perchloric acid. The vertical-rising sash may be raised to a maximum 28" operating height. The airflow requirements should be sized for the 28" operating height; if using sash stops then the airflow requirements can be reduced by approximately 40% at 18" or approximately 50% at 15".

S

Operating the Blower

Your Protector Fume Hood utilizes a remote style blower, which can be activated by turning the blower switch to "ON." You can validate the hood performance by watching smoke drawn into the hood face opening.

S

Operating the Lights

Your Protector Fume Hood utilizes a factory-wired fluorescent light to illuminate the hood interior. Simply turn the light switch to "ON" to operate.

S

Operating the Washdown Control Valve

Your Protector Radioisotope Fume Hood has a washdown control valve located on the upper left-hand side that controls water to the washdown spray nozzles. The nozzles are located behind the hood baffle and will washdown areas in the hood, which are inaccessible without removing the baffle. The washdown control system can also be set up to control the washdown rings used in your exhaust ductwork. This is done by connecting the water line to the plugged connection on

top of the washdown hood system or ordering additional washdown valves found in Chapter 2 and Chapter 7.

Working in Your Protector Radioisotope Fume Hood

Planning

- Thoroughly understand procedures and equipment required before beginning work.
- Arrange for minimal disruptions, such as room traffic or entry into the room while the hood is in use.

Start-up

- Turn on fluorescent light and hood blower.
- Slowly raise the sash.
- Check the baffle air slots for obstructions.
- Allow the hood to operate unobstructed for 5 minutes.
- Wear a long sleeved lab coat and rubber gloves. Use protective eyewear. Wear a protective mask if appropriate.

Loading Materials and Equipment

- Only load the materials required for the procedure. Do not overload the hood.
- High-level radioisotope materials (gamma particles) are not recommended for use in this Radioisotope Hood. Only use low-level alpha, lower level beta radioactive particles, and gaseous radioisotopes. Consult your Safety Officer should you have any questions.
- Do not obstruct the front air foil, or rear baffle slots.
- Large objects should not be placed close together and spaced above the liner bottom to permit airflow to sweep under the equipment.
- After loading the hood, wait one minute to purge airborne contaminants from the work area.

Work Techniques

- Keep all materials at least 6 inches inside of the sash, and perform all contaminated operations as far to the rear of the work area as possible. Segregate all clean and contaminated materials in the work area.
- Avoid using techniques or procedures that disrupt the airflow patterns of the hood.

Final Purging

- Upon completion of work, the hood should be allowed to operate for two to three minutes undisturbed, to purge airborne contaminants from the work area before shutting down blower.

Unloading Materials and Equipment

- Objects in contact with contaminated material should be surface decontaminated before removal from the hood.
- All open trays or containers should be covered before being removed from the hood.

Chapter 5: Using the Protector Radioisotope

Final Decontamination

- Water washdown or flooding the surfaces with water is standard practice and results in about 50% reduction of contamination on the 304 stainless steel liner material due to radiochemicals.
- A secondary means of surface decontamination involves the use of organic solvents that can result in up to 75% reduction of contamination due to radiochemicals.
- Even though 304 stainless steel seams are smooth to the touch, there are microscopic voids within the surface in which radioactive particles can lodge and are not easily removed except by abrasion decontamination. Abrasion decontamination consists of sanding and/or scouring operations to dislodge radioactive particles from the stainless steel hood liner.

Shutdown

- Only turn off the fluorescent light and hood blower if final washdown is complete, and then close the sash.

CHAPTER 6

MAINTAINING THE PROTECTOR

RADIOISOTOPE

Now that you have an understanding of how to work in the fume hood, we will review the suggested maintenance schedule and the common service operations necessary to maintain your fume hood for peak performance.

Service Safety Precautions

- If performing any electrical maintenance, always disconnect the power at the main disconnect.
- If performing decontamination inside the fume hood, consult your safety officer for proper personal protective equipment and procedure.
- Since some service operations require a step ladder, always use proper safety and consult your safety officer
- If performing maintenance on any service lines, always shut off the supply first.
- Some removable components may be heavy, follow safe-lifting guidelines.
- Verify all components are installed correctly with performance verified before conducting normal operations.

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Chapter 6: Maintaining the Protector Radioisotope

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
Routine Maintenance Schedule

Weekly

- Using ordinary dish soap to clean the surface inside of the fume hood, and the work surface.
- Using an appropriate glass cleaner, clean the sash and all glass surfaces.
- Operate the fume hood blower, noting the airflow velocity through the hood using a source of visible smoke.
- Decontaminate 304 stainless steel surfaces (if used daily) from radiochemicals as outlined in Chapter 5 Final Decontamination.

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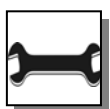Biannually

- The sash assembly should be checked to ensure proper operation and to make sure there are no signs of abnormal wear on the sash pulleys, cables and clamps.

Routine Service Operations

Front Panel Removal:

1. Simply lift the front panel up and then away from the hood to provide access to the top.

A**S**

Changing the Fluorescent Lamp:

1. Turn light switch to “OFF”.
2. Remove the front panel as noted earlier.
3. Reach over the front header of the hood and remove knock out plugs at both ends of fixture. Lift fixture up.
4. Remove the fluorescent lamp by pushing it out of the spring-loaded lamp socket and swinging it out of the other lamp socket.
5. Install the new lamp by reversing the removal procedure.

CHAPTER 7

MODIFYING THE PROTECTOR RADIOISOTOPE

There are several ways to modify the fume hood for your individual requirements. These include the addition of work tops, service fixtures, air monitor, electrical duplex outlets, ceiling enclosures, and rear panels.

Installing Hardboard Support Surface

Your Protector Radioisotope Fume Hood requires a hardboard support surface to properly support the integral work surface of the one-piece stainless steel liner. The hardboard support is shipped and included with the Radioisotope Fume Hood. Be sure to secure the hardboard support per the instructions in Chapter 3.

Installing Ceiling Enclosures Above the Fume Hood

Your Protector Fume Hood has mounting holes to accept a ceiling enclosure to close off the area between the top of the hood and the ceiling. Contact Labconco Customer Service for ordering information.

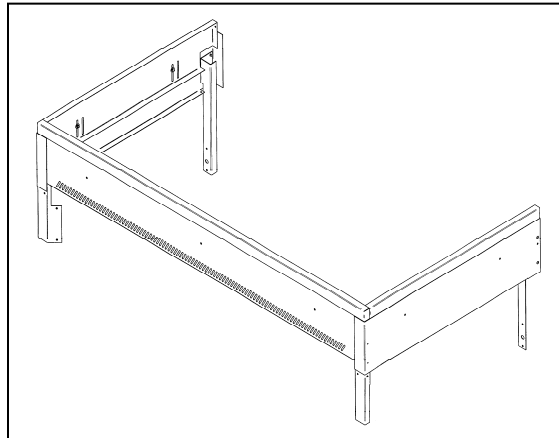
Figure 7-1

Installing Rear Panels Behind the Fume Hood

Your Protector Fume Hood can be modified to add a rear panel behind the fume hood when the fume hood is placed on an island (not available on Pass-Through Fume Hood). Contact Labconco Customer Service for ordering information.


Installing Guardian™ Digital Airflow Monitor or Guardian™ Airflow Monitor

The Guardian Airflow Monitors continuously monitor face velocity through the fume hood opening. The fume hood right corner post is factory prepared to mount either monitor. Contact Labconco Customer Service to order.

Sash Stop Kit – Field Installation (P/N 9410300)

The sash stop kit restricts how far a vertical-rising sash may be opened. This small plastic device may be easily field installed on the fixture corner post of any fume hood.

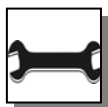Installing an Electrical Duplex Outlet

Your Protector Fume Hood can be ordered with duplex outlets, however, if you ordered a model without an electrical duplex outlet you can have one installed in the field by a qualified electrician. Contact Labconco Customer Service for ordering information. (Not acceptable on explosion-proof hoods).

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CHAPTER 8

TROUBLESHOOTING

Refer to the following table if your fume hood fails to operate properly. If the suggested corrective actions do not solve your problem, contact Labconco for additional assistance.

PROBLEM	CAUSE	CORRECTIVE ACTION
Remote blower and lights won't operate	Wires not connected at junction boxes or switches.	Check connection of switches.
		Check connection to control box on top of unit.
	Circuit breakers tripped in building electrical supply.	Reset circuit breakers.
Remote blower won't operate, but lights work	Blower wiring is disconnected.	Inspect blower wiring and switch.
	Belt broken.	Replace belt.
	Blower motor is defective.	Replace blower motor.
Fume hood blower operates but lights will not operate	Lamp not installed correctly.	Inspect lamp installation.
	Lamp is defective.	Replace lamp.
	Lamp circuit breaker in building is tripped.	Reset the lamp circuit breaker.

Chapter 8: Troubleshooting

PROBLEM	CAUSE	CORRECTIVE ACTION
Fume hood blower operates but lights will not operate	Lamp wiring is disconnected.	Inspect lamp wiring.
	Defective lamp ballasts.	Replace lamp ballasts.
Contaminants escape outside of fume hood	Improper user techniques for the fume hood.	See “Certifying the Hood” Chapter 3 and “Safety Precautions” Chapter 4 sections in the manual. (Ref. Appendix D).
	Restriction of the baffle air slots or – blockage of the exhaust outlet.	Remove baffles to ensure that all air slots and the exhaust outlet are unobstructed.
	External factors are disrupting the fume hood airflow patterns or acting as a source of contamination.	See “Location Requirements” Chapter 2, “Certifying the Hood” Chapter 3, and “Safety Precautions” Chapter 4 sections of this manual. (Ref. Appendix D)
	Fume hood has improper face velocity.	Have fume hood re-certified and check remote blower exhaust system. Hood should have average face velocity of 100 fpm.
Vertical sash no longer operates smoothly	Cable is frayed or plastic protection is damaged.	Inspect cable and replace cable if worn or damaged immediately; otherwise injury could result.
	Pulley bearing is damaged.	Replace pulley, bearing or add grease.
	Cable has slipped off the pulleys.	Re-install, cable must be replaced immediately if damaged.
	Weight has broken pulleys.	Replace weight pulleys.
Combination A-Style sash no longer operates smoothly	Horizontal-sliding glass panels have come off the tracks.	Re-install horizontal-sliding glass on tracks.

PROBLEM	CAUSE	CORRECTIVE ACTION
Electrical duplex outlets no longer have power	Vertical-rising sash frame is distorted.	Place horizontal sliding glass symmetrically and pull sash down to air foil. Straighten damaged frame.
	Cable is frayed or has slipped off the pulleys.	Re-install, cable must be replaced immediately if damaged.
	Wires not connected or faulty duplex.	Check wire connection or replace duplex.
	Circuit breakers tripped in building electrical supply.	Reset circuit breakers.
Service valves no longer operate	Faulty building supply.	Inspect building supply shut off valves and appropriate pressures below 40 PSI.
	Valve no longer operates.	Replace valve and check for leaks.
Spray Nozzles no longer work	Supply line or outlet line has leaks.	Inspect line for leaks and fix any leaking plumbing connections.
	Poor water supply	Inspect building water supply.
	Clogged nozzle	Replace spray nozzles.

APPENDIX A

PROTECTOR RADIOISOTOPE COMPONENTS


Illustration A-1 indicate the location of the following service parts:

Protector Radioisotope Replacement Parts

Item	Quantity	Part No.	Description
1A	1	9823700	Valve, Labconco (Water) 1/4" Compression Fitting
1B	1	9823701	Valve, Labconco (Water) 3/8" Compression Fitting
1C	1	9817000	Valve, Labconco 1/4" Compression Fitting (AIR, GAS, VAC, NIT, etc.)
1D	1	9817001	Valve, Labconco 3/8" Compression Fitting (AIR, GAS, VAC, NIT, etc.)
1E	1	9823702	Valve, Labconco Deionized 1/4" Compression Fitting
1F	1	9823703	Valve, Labconco Deionized 3/8" Compression Fitting
1G	1	9818000	Nut, Valve Mounting. (Labconco)
2A	1	9826800	WaterSaver Valve/Gooseneck -GRN
2B	1	9826801	WaterSaver Valve/Connector (VAC) – YEL
2C	1	9826802	WaterSaver Valve/Connector (AIR) – ORG
2D	1	9826803	WaterSaver Valve/Connector (GAS) – BLU
2E	1	9826805	WaterSaver Valve/Connector (HOT WATER) – RED
2F	1	9826806	WaterSaver Valve/Connector (CW) – GRN
2G	1	9826807	WaterSaver Valve/Connector (STEAM) – BLK
2H	1	9826808	WaterSaver Valve/Connector (NITROGEN) – BRN
2I	1	9826809	WaterSaver Valve/Connector (OXYGEN) – LIGHT GREEN
2J	1	9826810	Swivel Gooseneck only – GRN
2K	1	9826812	Swivel Gooseneck only – WHITE
3	1	9818700 thru 08	Knobs (GRAY, GRN, BLU, ORG, YEL, RED, WHT, BLK, BRN)
4A	1	9818800	Hose Barb, GRAY – (NEUTRAL OR ARGON) – NOT SHOWN
4B	1	9818801	Hose Barb, GREEN - (COLD WATER) – NOT SHOWN
4C	1	9818802	Hose Barb, BLUE – (GAS) – NOT SHOWN
4D	1	9818803	Hose Barb, ORANGE – (AIR) – NOT SHOWN
4E	1	9818804	Hose Barb, YELLOW – (VACUUM) – NOT SHOWN
4F	1	9818805	Hose Barb, RED – (HOT WATER) – NOT SHOWN
4G	1	9818806	Hose Barb, WHITE – (DEIONIZED WATER) – NOT SHOWN
4H	1	9818807	Hose Barb, BLACK – (NEUTRAL OR STEAM) – NOT SHOWN
4I	1	9818808	Hose Barb, BROWN – (NITROGEN) – NOT SHOWN
4J	1	9819000	Nut, Hose Barb – NOT SHOWN
5	1	9825500	Label, Knob (contains all the labels)
6	1	9818900	Lens, Knob
7A	1	9947100, 01, 02	115V Duplex Receptacle (GRAY) Right, Left 4' - 6', Left 8' w/ wires
7B	1	9818200	Cover Plate 115V Duplex
7C	1	9947103, 04, 05	115V GFCI Duplex Receptacle (GRAY) Right 4' - 6', Left 8' w/ wires
7D	1	9818100	Cover Plate, 115V GFCI
7E	1	9818300	Cover Plate, Blank
8A	1	9721901	Lamp, Fluorescent (T8 x 3') – use on 4' & 8' Hoods – NOT SHOWN
8B	1	9721900	Lamp, Fluorescent (T8 x 4') – use on 5' & 6' Hoods – NOT SHOWN

Appendix A: Protector Radioisotope Components

Item	Quantity	Part No.	Description
8C	1	9945300	Fluorescent Light Fixture, 36" – use on 4' Hoods
8D	1	9945301	Fluorescent Light Fixture, 48" – use on 5', 6' Hoods
8E	1	9945302	Fluorescent Light Fixture, 36" Combo – use on 8' Hoods
10A	1	1302301	Switch, Rocker
10B	1	1327500	Switch, Plug (Fills cutout when switch is not used)
11	1	9818400	Access Cover
12A	1	9409800	Side Panel, 24" internal deep hoods
12B	4	1885308	Screw, Machine #10-24 x .50 Phillips
13A	1	9466400	Eco-Foil 4'
13B	1	9466401	Eco-Foil 5'
13C	1	9466402	Eco-Foil 6'
13D	1	9466403	Eco-Foil 8'
14	4	1861400	Pulley, Front, 1-3/16 Dia.
15	2	4949902	Cable, Sash 130" – NOT SHOWN
16	2	9709300	Plastic Pulley, (Rear 2")
17	2	9713300	Bumper, Rubber – NOT SHOWN (lower sash bumper)
18	4	1934100	Bronze Bearing, Flanged Rear – NOT SHOWN
19	4	1920100	Clamp, Cable Replacement – NOT SHOWN
20	2	9935800	Threaded Connector – NOT SHOWN (to attach weight to cable)
21	1	9946300 or 01	Wiring Harness, Main
22	1	9409601, 02, 03, 05	Front Panel, 4', 5', 6', 8'

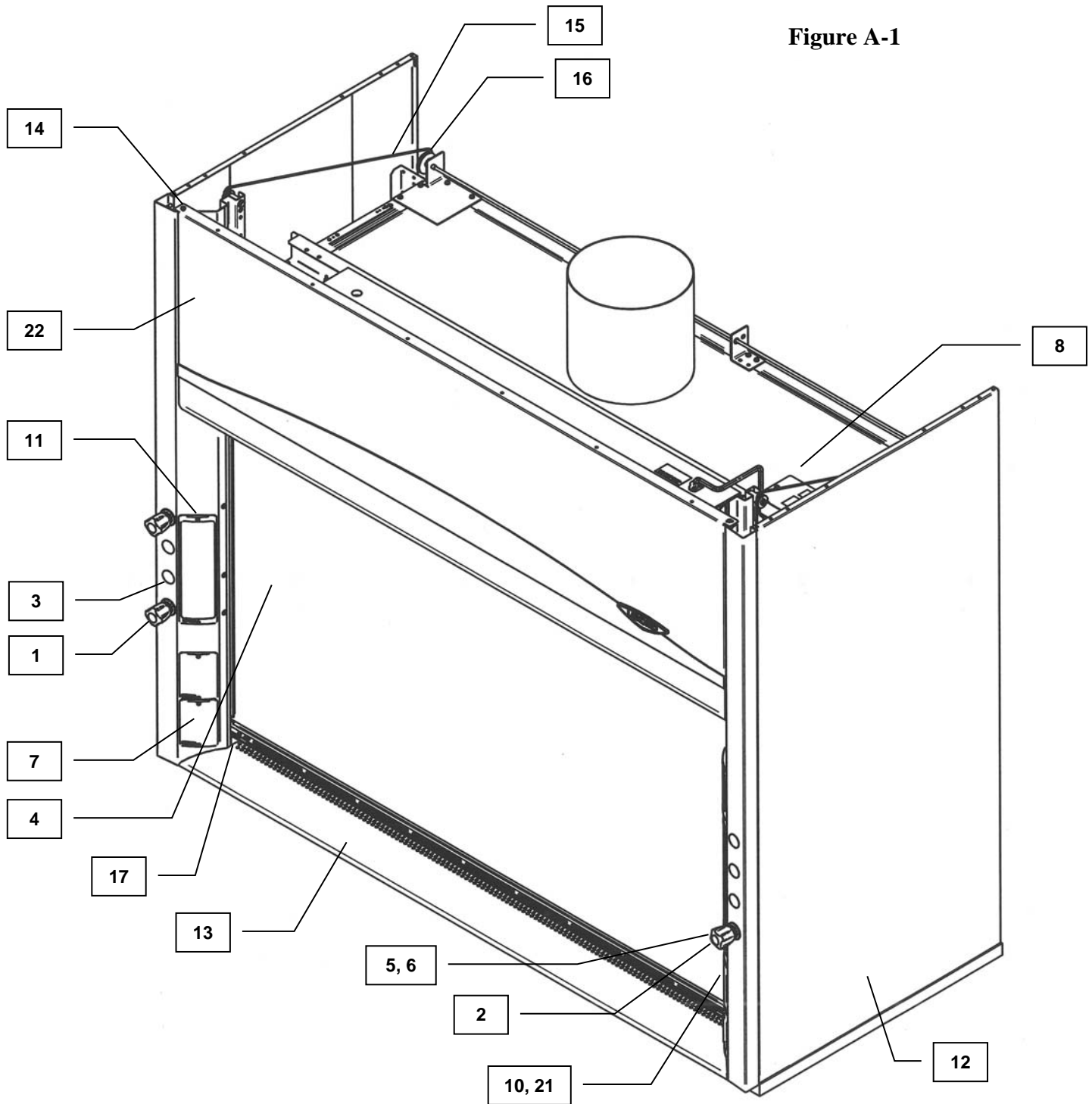Appendix A: Protector Radioisotope Components**Figure A-1**

APPENDIX B

PROTECTOR RADIOISOTOPE

DIMENSIONS

Radioisotope Model Widths

	4'	5'	6'	8'
A	48.00	60.00	72.00	96.00
B	38.25	50.25	62.25	86.25

Dimensions in inches. CAD layouts available upon request.


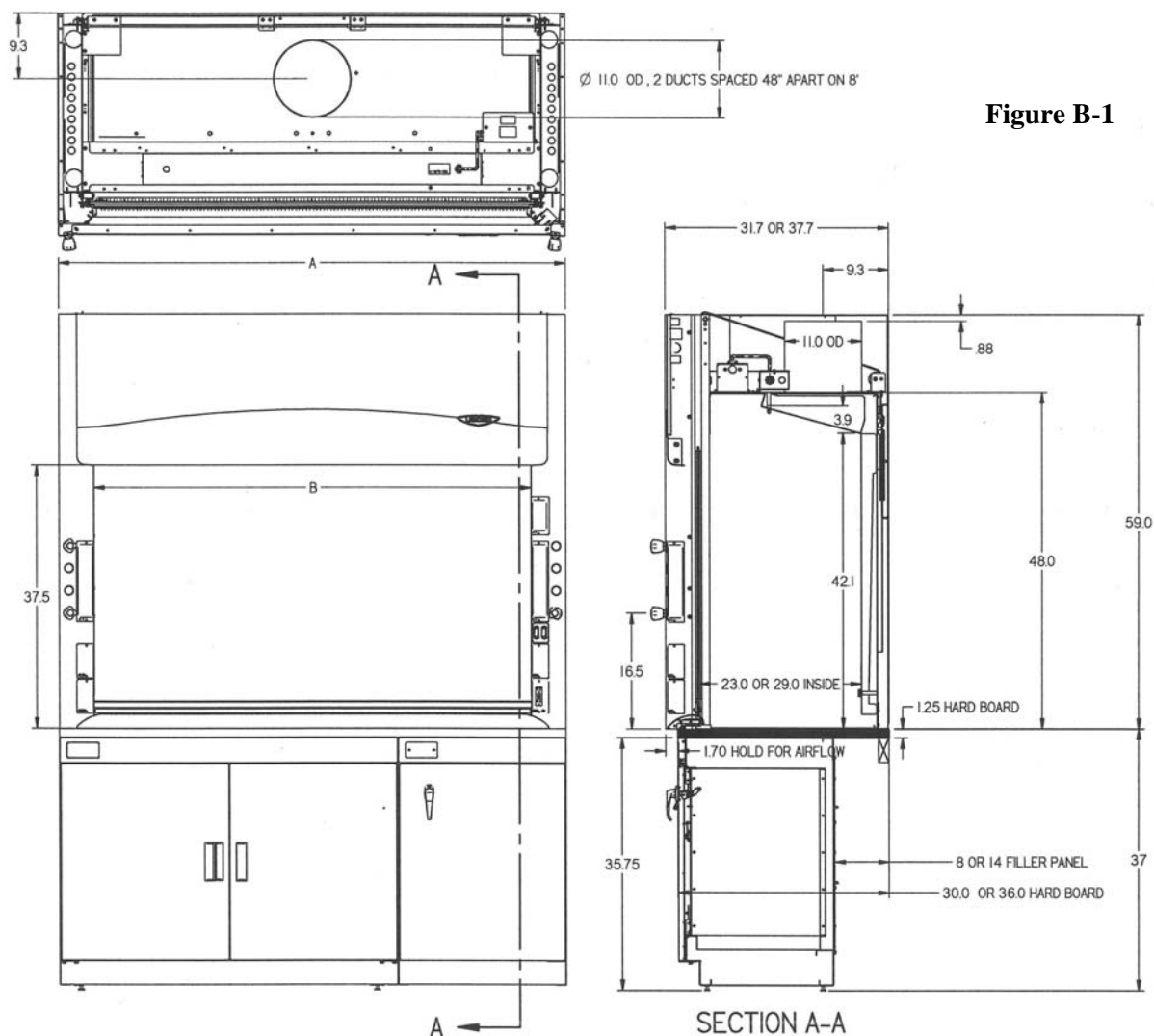
Figure B-1

Appendix B: Protector Radioisotope Dimensions**Figure B-2**

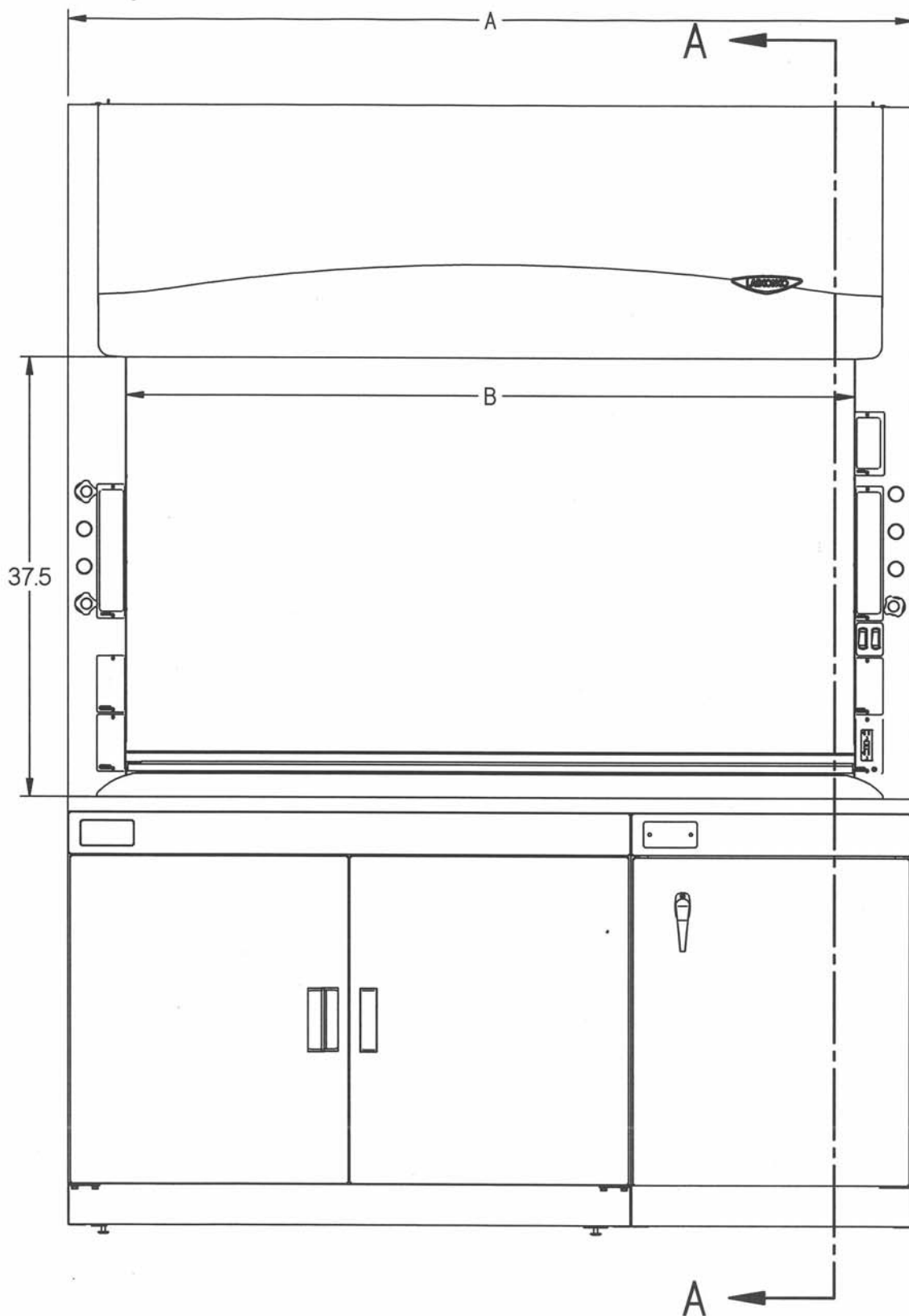Appendix B: Protector Radioisotope Dimensions

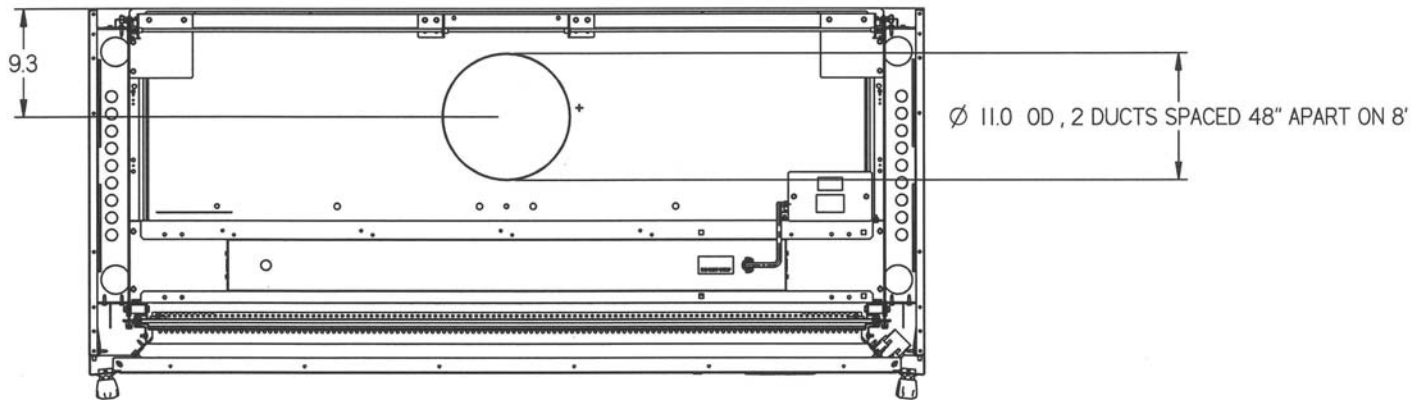
Figure B-4

APPENDIX C

PROTECTOR RADIOISOTOPE SPECIFICATIONS

Environmental Conditions


- Indoor use only.
- Maximum altitude: 10,000 feet (3,048 meters).
- Ambient temperature range: 41° to 104°F (5° to 40°C).
- Maximum relative humidity: 80% for temperatures up to 88°F (31°C), decreasing linearly to 50% relative humidity at 104°F (40°C).
- Main supply voltage fluctuations not to exceed $\pm 10\%$ of the nominal voltage.
- Transient over-voltages according to Installation Categories II (Over-voltage Categories per IEC 1010). Temporary voltage spikes on the AC input line that may be as high as 1500V for 115V models and 2500V for 230V models are allowed.
- Used in an environment of Pollution degrees 2 (i.e., where normally only non-conductive atmospheres are present). Occasionally, however, a temporary conductivity caused by condensation must be expected, in accordance with IEC 664.
- Electrical Ratings

<u>Volts AC</u>	<u>Phase</u>	<u>Cycle</u>	<u>AMP*</u>
115	1	50/60	10
115	1	50	10
115	1	60	10
230	1	50/60	5
230	1	50	5
230	1	60	5

*Does not include current rating of receptacles.

Appendix C: Protector Radioisotope Specifications

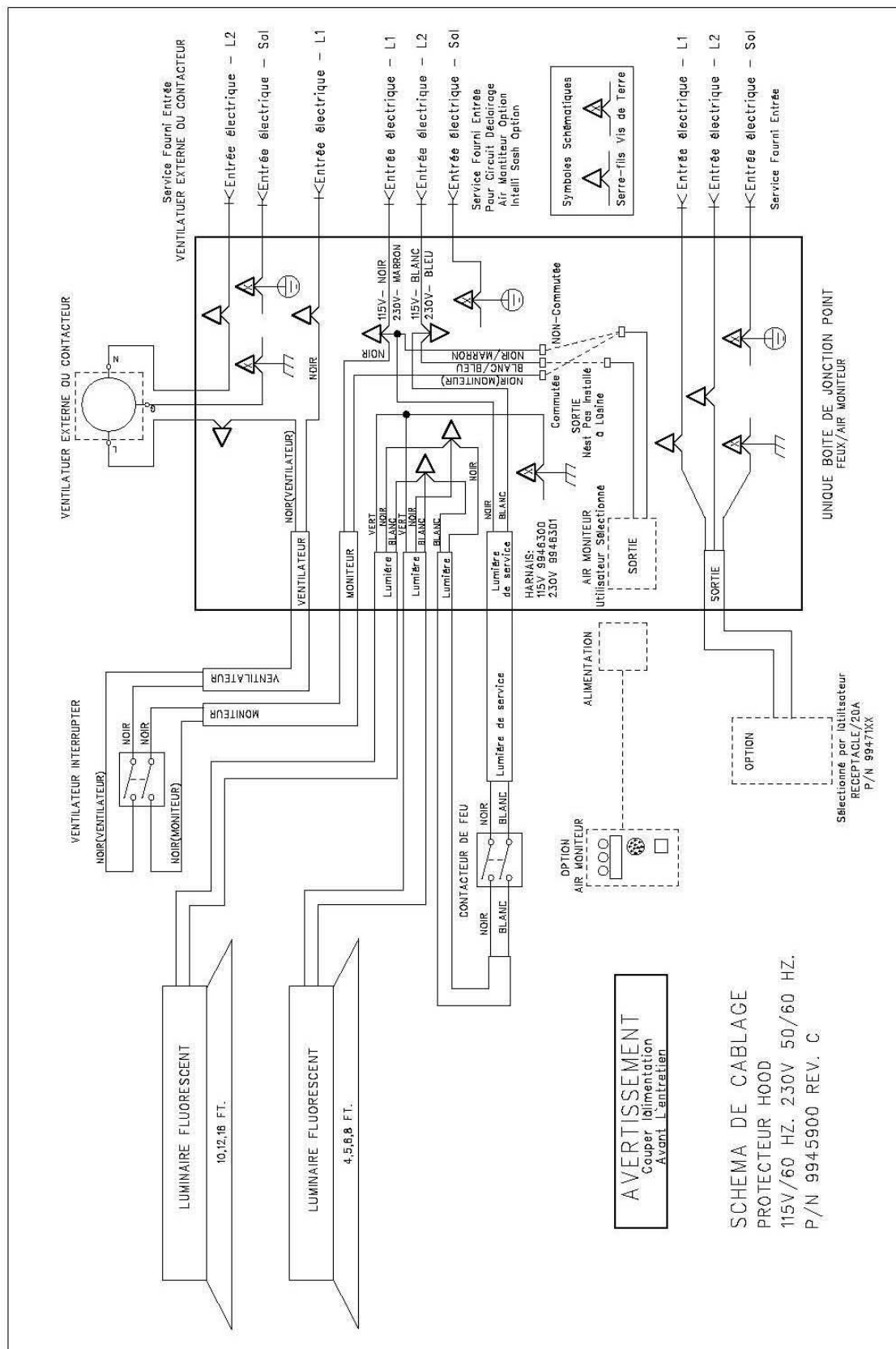APPENDIX D

SERIAL NUMBER TAG

DESCRIPTION

Serial tag includes standard information with the following changes:

The portion of the tag for electrical information is labeled with the following units: Volts AC, Phase, Cycle, AMP.

Each hood will include one of the following:

<u>Volts AC</u>	<u>Phase</u>	<u>Cycle</u>	<u>AMP</u> 
115	1	50/60	10X
115	1	50	10X
115	1	60	10X
230	1	50/60	5X
230	1	50	5X
230	1	60	5X

Where  may be a separate label indicating: Caution, See Manual.

The X above will be one of several possible alpha characters as defined in the manual as follows:

- X: No additional duplexes on hood, noted amperage applies.
- B: This fume hood includes **one** additional electrical receptacle, individually wired to the field wired box and individually rated as 115V, single phase, 60Hz, 20 Amps. Each duplex can be wired on a dedicated circuit rated at 20A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 20A electrical service they are wired to.
- C: This fume hood includes **two** additional electrical receptacles, individually wired to the field wired box and individually rated as 115V, single phase, 60Hz, 20 Amps. Each duplex can be wired on a dedicated circuit rated at 20A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 20A electrical service they are wired to.

- D: This fume hood includes **three** additional electrical receptacles, individually wired to the field wired box and individually rated as 115V, single phase, 60Hz, 20 Amps. Each duplex can be wired on a dedicated circuit rated at 20A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 20A electrical service they are wired to.
- E: This fume hood includes **four** additional electrical receptacles, individually wired to the field wired box and individually rated as 115V, single phase, 60Hz, 20 Amps. Each duplex can be wired on a dedicated circuit rated at 20A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 20A electrical service they are wired to.
- F: This fume hood includes **one** additional electrical receptacle, individually wired to the field wired box and individually rated as 230V, single phase, **60Hz**, 10 Amps. Each duplex can be wired on a dedicated circuit rated at 10A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 10A electrical service they are wired to.
- G: This fume hood includes **two** additional electrical receptacles, individually wired to the field wired box and individually rated as 230V, single phase, **60Hz**, 10 Amps. Each duplex can be wired on a dedicated circuit rated at 10A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 10A electrical service they are wired to.
- H: This fume hood includes **three** additional electrical receptacles, individually wired to the field wired box and individually rated as 230V, single phase, **60Hz**, 10 Amps. Each duplex can be wired on a dedicated circuit rated at 10A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 10A electrical service they are wired to.
- J: This fume hood includes **four** additional electrical receptacles, individually wired to the field wired box and individually rated as 230V, single phase, **60Hz**, 10 Amps. Each duplex can be wired on a dedicated circuit rated at 10A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 10A electrical service they are wired to.
- K: This fume hood includes **one** additional electrical receptacle, individually wired to the field wired box and individually rated as 230V, single phase, **50Hz**, 10 Amps. Each duplex can be wired on a dedicated circuit rated at 10A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 10A electrical service they are wired to.
- L: This fume hood includes **two** additional electrical receptacles, individually wired to the field wired box and individually rated as 230V, single phase, **50Hz**, 10 Amps. Each duplex can be wired on a dedicated circuit rated at 10A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 10A electrical service they are wired to.

Appendix D: Serial Number Tag Description

- M: This fume hood includes **three** additional electrical receptacles, individually wired to the field wired box and individually rated as 230V, single phase, **50Hz**, 10 Amps. Each duplex can be wired on a dedicated circuit rated at 10A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 10A electrical service they are wired to.
- N: This fume hood includes **four** additional electrical receptacles, individually wired to the field wired box and individually rated as 230V, single phase, **50Hz**, 10 Amps. Each duplex can be wired on a dedicated circuit rated at 10A or the duplexes may be ganged together on the same circuit with the total load of the ganged duplexes not exceeding the 10A electrical service they are wired to.

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APPENDIX E

REFERENCES

Many excellent reference texts and booklets are currently available. The following is a brief listing:

Laboratory Ventilation Standards

Federal Register 29 CFR Part 1910

Non-mandatory recommendations from "Prudent Practices".

- Fume hoods should have a continuous monitoring device
- Face velocities should be between 60-100 linear feet per minute (lfpm)
- Average 2.5 linear feet of hood space per person

Occupational Health and Safety

U.S. Department of Labor

200 Constitution Avenue N.W.

Washington, DC 20210

(202) 523-1452

Industrial Ventilation-ACGIH

- Fume hood face velocities between 60-100 lfpm
- Maximum of 125 lfpm for radioisotope hoods
- Duct velocities of 1000-2000 fpm for vapors, gasses and smoke
- Stack discharge height 1.3-2.0 x building height
- Well designed fume hood containment loss, <0.10 ppm

Industrial Ventilation, A Manual of Recommended Practice.

24th Edition, 2001

American Conference of Governmental Industrial Hygienists

1330 Kemper Meadow drive

Cincinnati, OH 45240-1634

(513) 742-2020

ASHRAE 110-1995 Method of Testing Performance of Fume Hoods

Evaluates fume hood's containment characteristics

- Three part test: Smoke generation, Face velocity profile, Tracer gas release @ 4 liters per minute
- Rated As Manufactured (AM), As Installed (AI) and As Used (AU)

American Society of Heating, Refrigerating, and Air Conditioning Engineers

1791 Tullie Circle N.E.

Atlanta, GA 30329

(404) 636-8400

ANSI Z9.5-2011 Laboratory Standard

Covers entire laboratory ventilation system.

- New and remodeled hoods shall have a monitoring device
- Ductless hoods should only be used with non-hazardous materials

American Industrial Hygiene Association

2700 Prosperity Avenue, Suite 250

Fairfax, VA 22031

(703) 849-8888

SEFA 1-2002

- Fume hood face velocities based on toxicity levels of chemicals
 - Class A – 125 to 150 fpm
 - Class B – 80 to 100 fpm
 - Class C – 75 to 80 fpm
- Test method – face velocity profile and smoke generation

Scientific Equipment & Furniture Association

1028 Duchess Drive

McLean, VA 22102

(703) 538-6007

NFPA 45 – 2002 Fire Protection for Laboratories Using Chemicals

- Laboratory hoods should not be relied on for explosion protection
- Exhaust air from fume hoods should not be recirculated
- Services should be external to the hood
- Canopy hoods only for non-hazardous applications
- Materials of construction should have flame spread of 25 or less
- 80 to 120 fpm to prevent escape

NFPA 30 – 2000 Flammable and Combustible Liquids Code

- Approved cabinets may be metal or wood
- Vent location on cabinets are required
- Venting of cabinets not a requirement

National Fire Protection Association

1 Batterymarch Park

P.O. Box 9101

Quincy, MA 02269-9101

(800) 344-3555

General References

American Conference of Governmental Industrial Hygienists. *Industrial Ventilation, A Manual of Recommended Practice*, Cincinnati, OH

ASHRAE Standard Committee. *ASHRAE Standard* Atlanta: ASHRAE Publications Sales Department, 1995

British Standards Institution, *Laboratory Fume Cupboards*. Parts 1, 2 and 3, London: 1990

Appendix E: References

Department of Labor, Occupational Safety and Health Administration, *29 CFR Part 1910, Occupational Exposures to Hazardous Chemicals in Laboratories, Final Rule*. Vol. 55, No. 21. Washington D.C.:1990

DiBerardinis, L. et al. *Guides for Laboratory Design, Health and Safety Considerations*. Wiley & Sons, 1987

McDermott, Henry, *Handbook of Ventilation for Contaminant Control, 2nd Edition*. Butterworth Publishers, 1985.

Miller, Brinton M. et al. *Laboratory Safety: Principles and Practices*. American Society for Microbiology, Washington, D.C.: 1986

NIH Guidelines for the Laboratory Use of Chemical Carcinogens. NIH Publication No. 81-2385.

Rayburn, Stephen R. *The Foundation of Laboratory Safety, A Guide for the Biomedical Laboratory*. Springer-Verlag, New York: 1990

Sax, N. Irving and Lewis, JR., Richard J. *Rapid Guide to Hazardous Chemicals in the Workplace*. Van Nostrand Reinhold, 1987.

Schilt, Alfred A. *Perchloric Acid and Perchlorates*. The G. Frederick Smith Chemical Company, Columbus, OH: 1979.

Steere, Norman. *CRC Handbook of Laboratory Safety, 2nd Edition*. CRC Press, 1971.

Protecting your
laboratory environment

LABCONCO

LABCONCO CORPORATION
8811 Prospect Avenue
Kansas City, MO 64132-2696

(816) 333-8811 • (800) 821-5525
FAX: (816) 363-0130
www.labconco.com

CE Declaration of Conformity

We: Labconco Corporation
8811 Prospect Avenue
Kansas City, MO 64132 USA
Ph. 816-333-8811 Fax. 816-363-0130
labconco@labconco.com

being the representative party in the EU solely responsible for this certification, hereby declare that the following product:

Equipment: Laboratory Equipment – Protector® Stainless Steel Laboratory Fume Hoods

Type, Model:	1204 Series	4' Protector Stainless Steel Radioisotope Hood
	1205 Series	5' Protector Stainless Steel Radioisotope Hood
	1206 Series	6' Protector Stainless Steel Radioisotope Hood
	1208 Series	8' Protector Stainless Steel Radioisotope Hood

Manufacturer: Labconco Corporation
8811 Prospect Avenue
Kansas City, MO 64132 USA
Ph. 816-333-8811 Fax. 816-363-0130
labconco@labconco.com

is in conformity with the standards listed below:

EN61010-1
EN61326-1
EN55022
EN61000-3-2
EN61000-3-3

following the provisions of the following directives:

2006/42/EC
2004/108/EC
2006/95/EC

when installed and operated in accordance with the manufacturers installation and operating instructions.

 KANSAS CITY USA 14 Nov 2013
Mark Schmitz Place and date of issue
Vice President, Research and Engineering
Labconco Corporation

Labconco P/N 3696095, Rev. -, ECO #H847